


Aa
Bb
Cc
Dd
Ee
Ff
Gg
Hh
Ii
Jj
Kk
Ll
Mm
Nn
Oo
Pp
Qq
Rr
Ss
Tt
Uu
Vv
Ww
Xx
Yy
Zz

Draining Rack, Wall Mounting

Single panel drying racks for wall mounting. Manufactured in white polythene coated mild steel.

Catalogue No	Type	Dimensions (hwx), mm
D602-321	Single wall mounted	510 x 480
D602-325	Single wall mounted	400 x 300
D602-329	Double panel free standing	510 x 480

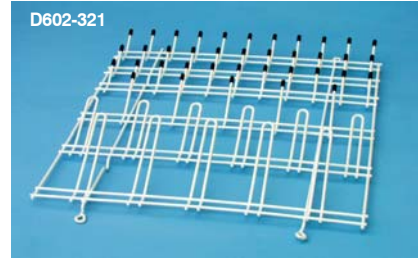

Draining Rack, High Impact Polystyrene

The 72 place support is attached to a back plate made from a single mould of 4mm rigid HIPS, and is stain and rust free. The peg sockets have closed ends eliminating leakage and consequent biohazards. The seamless drip channel collects waste which drains via drain tube provided. The 72 pegs are easily removable for cleaning, as well as to accommodate items of varying shapes. This item comes complete with mounting kit, draining hose and fixing template. Dimensions, 450 x 630mm.

Catalogue No
D602-340



Accessories: Draining Rack, Pegs

Pegs that are 6mm in diameter and can be attached to draining racks for small items such as test tubes. Supplied in packs of 11.

Catalogue No
D602-342



Drying Cabinet



Chalice drying cabinets are ideal for drying glassware and instruments or for applications requiring controlled warming. Our range has the following features:

- Constructed of single skin mild steel and finished in easy to clean stove powder paint for durability or stainless steel
- Adjustable, chromium plated grid wire shelves
- Heated through incoloy sheathed elements positioned below the cabinet floor
- Doors are made of toughened glass
- 113 and 226 litre models are fitted with sliding doors
- 540 litre model is fitted with hinged doors and is available with fan assisted circulation and/or variable rate extraction

Catalogue No	D693-513	D693-517	D693-521
Mild Steel			
Catalogue No	D693-514	D693-518	-
Stainless Steel			
Volume approx. litres	113	226	540
Ext. dims. (hwxwd), mm	559x730x360	666x928x449	1664x583x610
Int. Dims. (hwxwd), mm	460x728x341	567x926x430	1550x581x600
Maximum temperature, °C	50°C		
No of shelves	3	3	3
Shelf size (wxd), mm	810x330	910x415	565x585
Power rating, Watts	750	1000	1750
Net weight, Kg	30	40	80



Optional extras available on our 540 litre model

Option	Description
D693-537	Fan assisted circulation for improved drying performance (increases weight by 3Kg)
D693-538	Variable rate extractor for exhausting moisture (increases depth by 75mm and weight by 4Kg)
D693-539	Increase in maximum temperature to 65°C above ambient

